

Part 1 - 100# White Tag - 4" x 13 1/2" ($\pm 1/16"$) - 102mm x 342.9mm
Prints Black and Purple 7680U Inks
Laser 3 of 9 Barcode

VO Perf
Does Not Print

Glue Line Between Parts 1 & 4 Indicated in Gray
Does Not Print - Use Non-tinted Std glue


Fold Perf
Does Not Print

Part 2 - 903 Lot W231
4" x 2 3/32" (±1/16") - 102mm x 53.181mm
Prints Bio Black 586 Ink
Circles are 11mm ID

CUSTOMER	
APPROVED ☐	
NOT APPROVED ☐	
SIGNATURE	
NAME:	
DATE:	
EBF	
REF. 10550608	REV. AA
DATE:	
SIGNATURE	

Note: This PDF form layout is produced to a 1:1 scale. All copy and construction features are shown in their proper position per your specifications. Production variances will result in a potential $\pm 1/16"$ (1.6mm) tolerance.

3 Glue lines between Parts 2 & 3 indicated in Green. Do Not Print

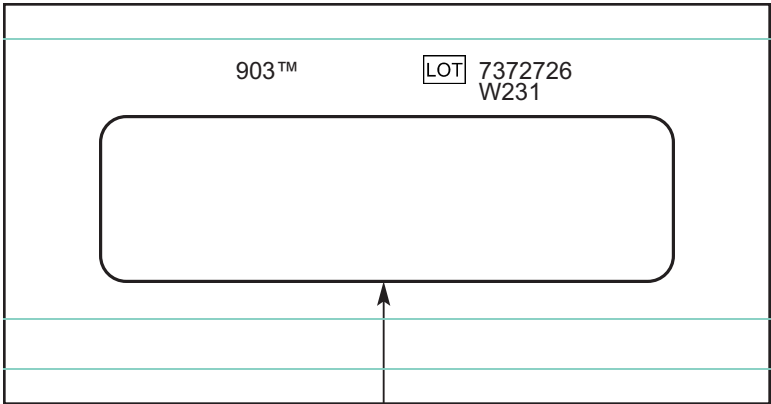

Die cut is 3" x 7/8" and
does not print

Part 3 Backer - 100# White Tag
4" x 2 3/32" (±1/16") - 102mm x 53.181mm
Prints Black Ink

903 DBS Cassette Card - Alpha Detect
10550608 Rev. AA
Job # 7372726-011
07-29-26
Fourth Proof

CUSTOMER	
APPROVED ☐	
NOT APPROVED ☐	
SIGNATURE	
NAME:	
DATE:	
EBF	
REF. 10550608	REV. AA
DATE:	
SIGNATURE	

Note: This PDF form layout is produced to a 1:1 scale. All copy and construction features are shown in their proper position per your specifications. Production variances will result in a potential ± 1/16" (1.6mm) tolerance.

Part 4 - White Self Contained - 4" x 13 1/2" (±1/16") - 102mm x 342.9mm - Prints Black Ink

Fold Perf
Does Not Print

Part 5 - 15# White Bond
4" x 2" ($\pm 1/16"$) - 102mm x 50.8mm
No Print

CUSTOMER	
APPROVED ☐	
NOT APPROVED ☐	
SIGNATURE	
NAME:	
DATE:	
EBF	
REF. 10550608	REV. AA
DATE:	
SIGNATURE	

Note: This PDF form layout is produced to a 1:1 scale. All copy and construction features are shown in their proper position per your specifications. Production variances will result in a potential $\pm 1/16"$ (1.6mm) tolerance.

INFORMATION ON FORM must be completed for identification and reporting of results.
Use Ballpoint Pen ONLY.

For questions please see lab contact information on outer box.

BLOOD COLLECTION INSTRUCTIONS

- Wash hands with warm water and dry thoroughly
- With Alcohol Prep, clean the tip of the finger and allow to AIR DRY, approximately 12-15 seconds. DO NOT BLOW ON THE FINGER
- Remove the protective blue tip from the lancet and dispose of protective tab immediately in trash.
- Warning: choking hazard. Keep away from children, device contains small parts. For more information on lancet use, visit www.bd.com.
- Hold hand palm up. Place tip of lancet slightly off center of the finger tip. Quickly and firmly press the lancet until it "clicks", then release
- When blood begins to flow, wipe the first drop away with a sterile gauze pad. If blood does not readily flow, apply pressure to the base of the finger and massage toward the tip. The drops should be large and "hanging" before being touched to the center of each circle.
- DO NOT PRESS THE FINGER TO THE CARD
- If blood flow begins to diminish, follow the initial blood collection procedures in order to fill the remaining circles
- Fill each circle with a single large drop of blood. Blood should completely cover the inside of each circle. If a drop is too small to completely fill a circle, additional small drops may be placed next to each other, BUT TRY NOT TO OVERLAP THE DROPS
- After all the circles have been filled, use a sterile gauze pad to stop the blood flow and then apply the bandage

IT IS VERY IMPORTANT TO ALLOW THE SAMPLES TO COMPLETELY AIR DRY AT ROOM TEMPERATURE BEFORE BEING SHIPPED.
FOR PROPER DRYING, STAND CARD ON ITS SIDE FOR AT LEAST 3 HOURS (OR OVERNIGHT). DO NOT REFRIGERATE. DO NOT
STACK WET CARDS ON TOP OF ONE ANOTHER TO PREVENT CONTAMINATION OF THE SAMPLES.

2031-07-31



10550608
REV. AA

VO Perf
Does Not Print

Fold Perf
Does Not Print

903 DBS Cassette Card - Alpha Detect
10550608 Rev. AA
Job # 7372726-011
07-29-26
Fourth Proof

CUSTOMER	
APPROVED ☐	
NOT APPROVED ☐	
SIGNATURE	
NAME:	
DATE:	
EBF	
REF. 10550608	REV. AA
DATE:	
SIGNATURE	

Note: This PDF form layout is produced to a 1:1 scale. All copy and construction features are shown in their proper position per your specifications. Production variances will result in a potential ± 1/16” (1.6mm) tolerance.